aquarick

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- ④ The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
- ⑤ If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
- ⑥ Never operate this product if it has a damaged cord, plug or charger base. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue using if the product appears damaged in any way.
- ⑦ Keep the toothbrush, charger base and power cord away from heated surfaces.
- ⑧ If you experience discomfort or pain in your mouth, discontinue the use of this product and contact a dentist.
- ⑨ Do not attempt to remove the battery.
- ⑩ Do not clean the toothbrush in the dishwasher.
- ⑪ Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

4. MEDICAL-WARNINGS


- ① If you have had oral surgery or gum surgery within the past two months, consult your dental provider before using the product.
- ② Consult your dental provider if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- ③ Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical devices.
- ④ Consult your physician prior to using the toothbrush if you have any medical concerns.
- ⑤ The toothbrush has been tested and is certified in accordance with CE standards.
- ⑥ Never use while sleeping or drowsy.

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
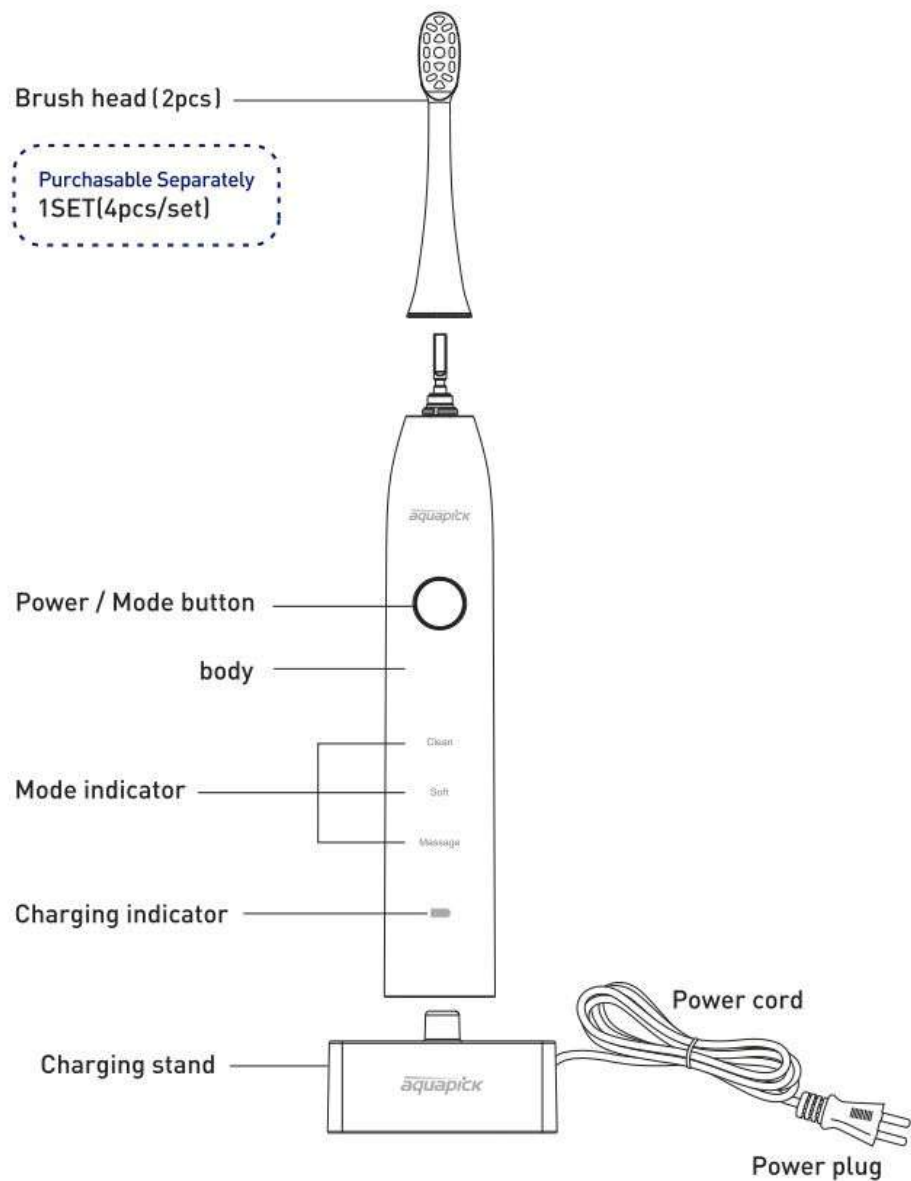
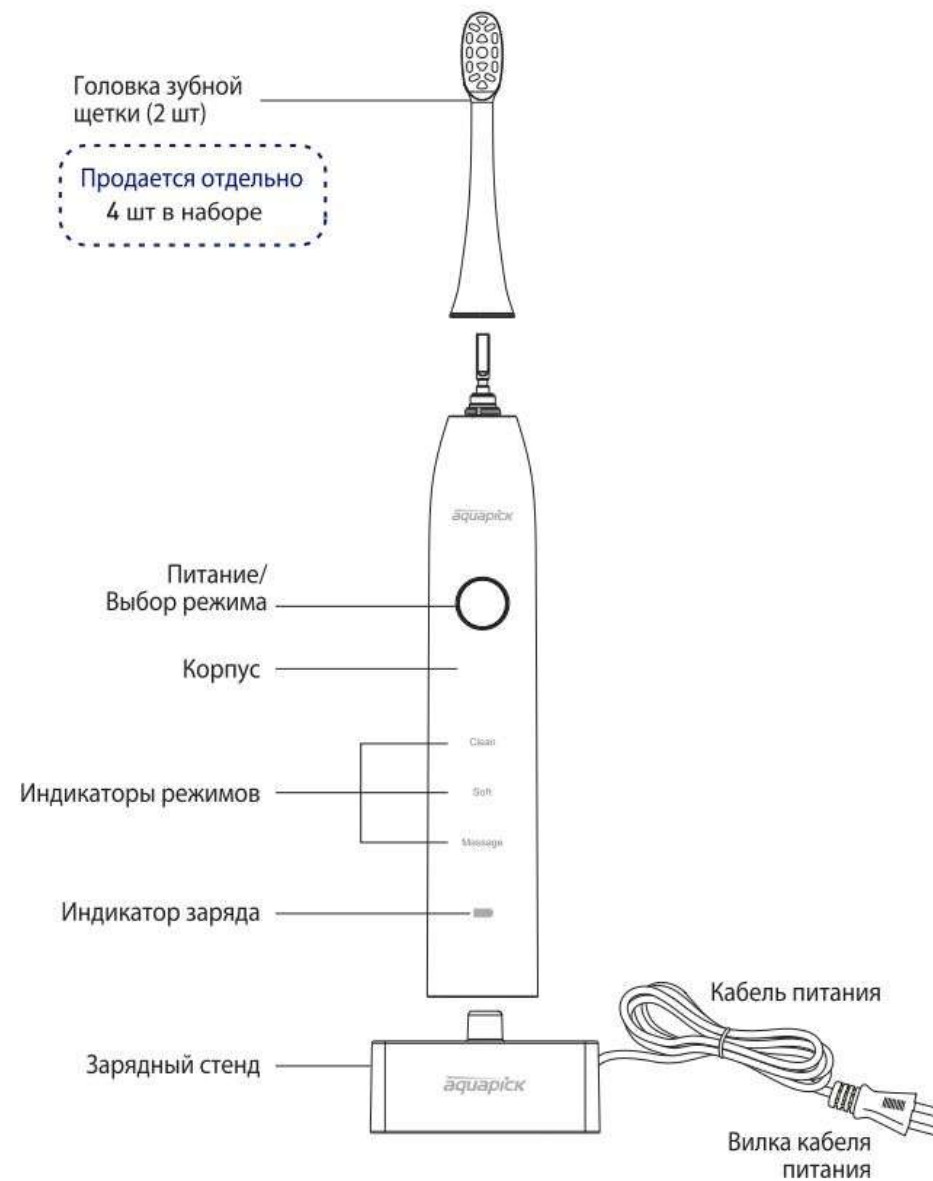
GETTING STARTED

Charging the toothbrush

Be sure to turn off the toothbrush before charging. Place the handle in the charger. The charge indicator blinks when the device is charging.

Notes:

- 1 It takes 24 hours to charge the battery fully.
- 2 The toothbrush must be charged for 24 hours, if:
 - a. this is the first time to use the toothbrush;
 - b. it has not been used for more than three months;
- 3 To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.
- 4 During charging, the toothbrush will get warm, but this is normal.
- 5 The operating duration will become shorter as the battery gets older. If the toothbrush operates only for short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact service center for support.

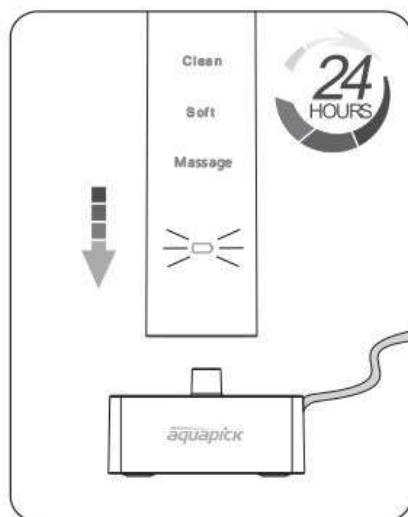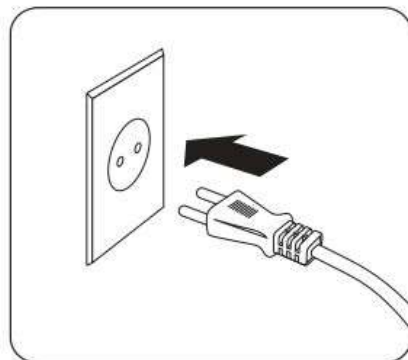The remaining charge indicator

When charging

- 1 Below 10%: Red lamp flickers
- 2 Over 10%: Green lamp flickers
- 3 Full charge: Green lamp always on.

When operating

- 1 Below 10%: Red lamp always on
- 2 Below 10%: Green always on

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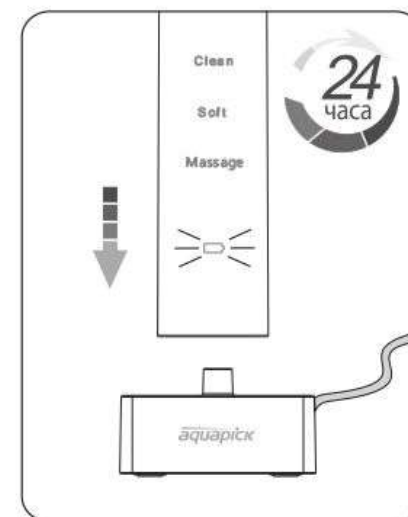
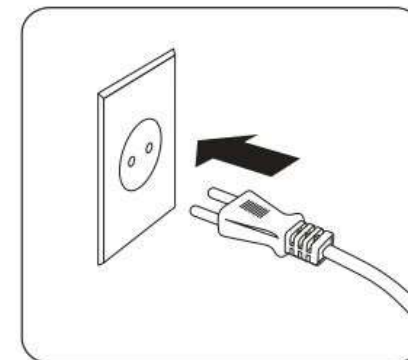
HOW TO USE

- Slide the brush head into the handle. Push brush head down until it locks into place on handle.
Note: There is small gap between silver ring and the handle. (Fig 1)
- Wet the bristles and apply a small amount of toothpaste.
- Place the bristles against the teeth at a slight angle towards the gum line. (Fig 2)
- Press the Power/Mode button to turn on the toothbrush and it starts working at Clean mode.




The brushing modes are as follows:

- Clean Mode:** Standard mode for superior teeth cleaning
- Soft Mode:** Gently cleans your teeth or any other sensitive areas;
- Massage Mode:** Pulsates for gentle gum stimulation.

During brushing, if the Power/Mode button is pressed, the toothbrush will work at next mode and start another 2 minutes cycle and the mode indicator will light up.

After pressing the Power/Mode button at first, if the button is pressed within 5 seconds again, its mode is changed and if the button is pressed after 5 seconds, the power is turned off.

- Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.
- In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. (Fig.3)

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
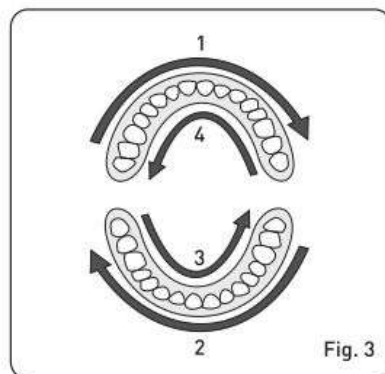
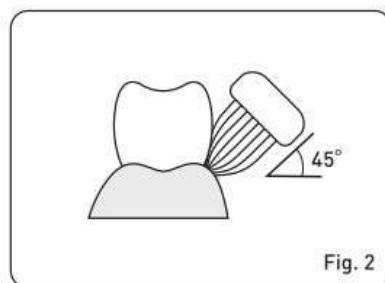
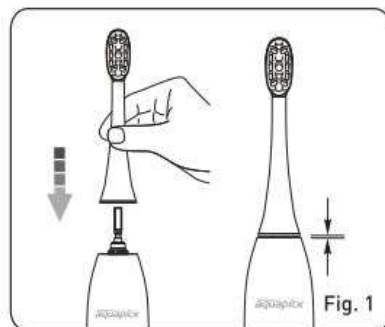
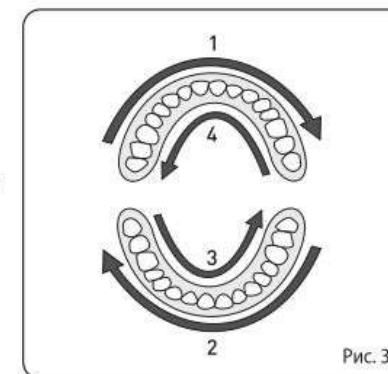
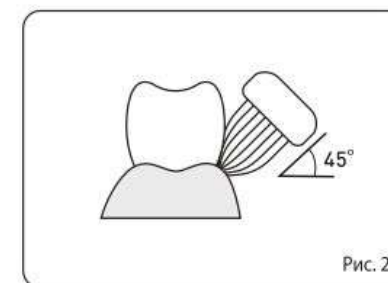
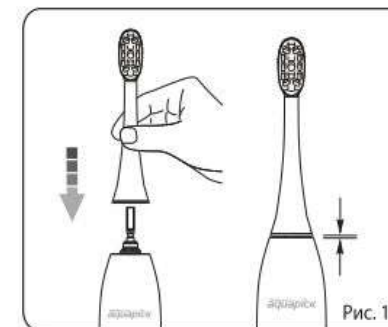
The brushing action will pause between each 30 second interval, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use.

- 7 Pull brush head out from the handle. (Fig 4)
- 8 Rinse both the brush head and handle in running water and dry them. (Fig 5)
- 9 Put the toothbrush back into the charger.

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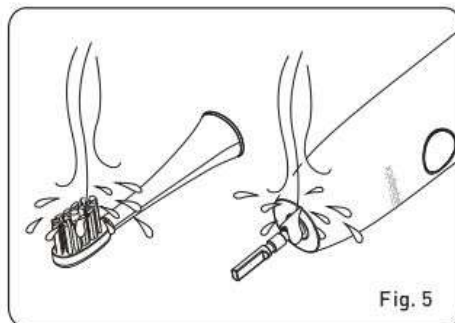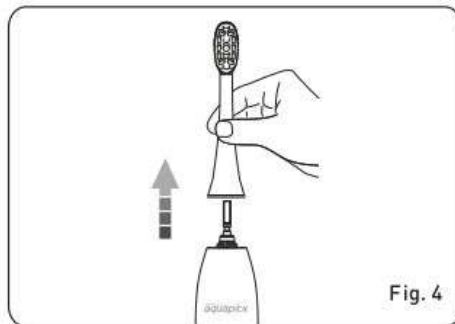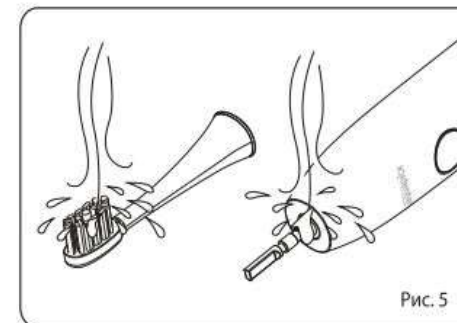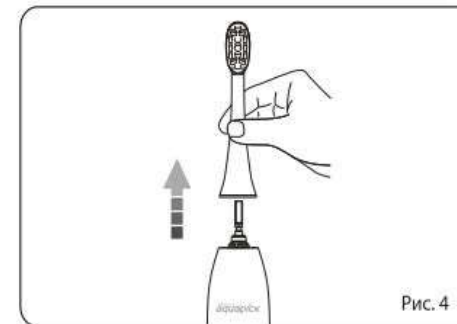MAINTENANCE

Toothbrush

- 1 Remove the brush head.
- 2 Wipe all surface areas with a damp cloth and the toothbrush back into the charger.
※ Do not submerge the toothbrush into water.
- 3 Rinse the brush head and bristles after each use.
- 4 Do not clean in the dishwasher.

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

Charger

- 1 Using a soft cloth, wipe off any dirt on the charger and the storage.
- 2 Clean the plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheat or fire.

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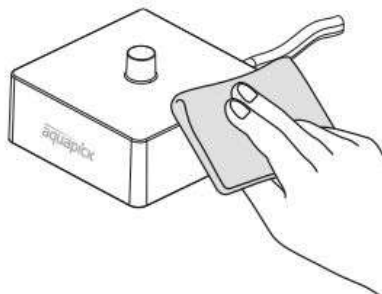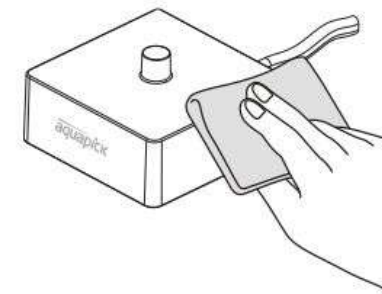TROUBLE SHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 24 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 24 hours.

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DISPOSAL

Management requirements.




If you have questions about recycling, contact your local waste management office.

Removing and recycling the battery

Do not disassemble the product arbitrarily except for product disposal.

Doing so might destroy the product. Observe basic safety precautions.

- 1 Before removing the battery, be sure to completely drain the unit of power.
- 2 Use a screwdriver to remove the screw cover and battery door.
- 3 Hold the handle upside down, press the toothbrush shaft, and remove the two hooks to release the internal components of the handle.
- 4 Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.
- 5 Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

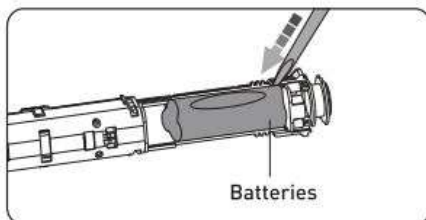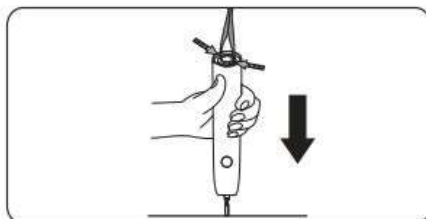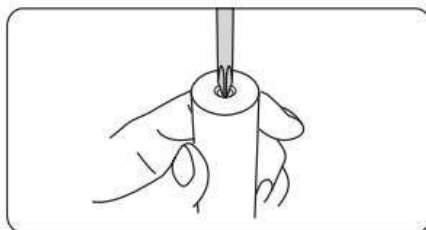Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to mains again after you have removed the battery.

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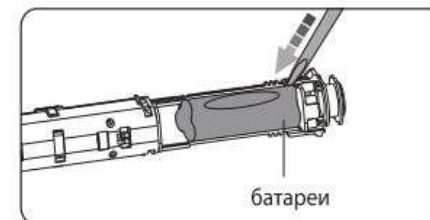
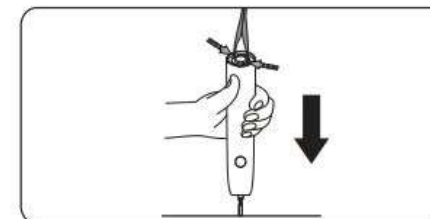
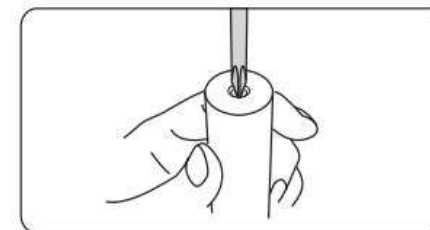
IMPORTANT SAFE GUARDS

The Sonic electric toothbrush is safe and effective when used according to directions. Like any electrical appliance, however, basic safety precautions should be followed, especially when children are present.

READ ALL INSTRUCTIONS BEFORE USING:

1. DANGER-WATER

To prevent the risk of burns, electrocution, fire, or injury to persons:

- 1 Do not place or store the charger base where it can fall or be pulled into a bathtub, sink or toilet.
- 2 Do not place or drop the charger into water or other liquid.
- 3 Do not use while bathing.
- 4 Unplug charger base before cleaning and assure that it is dry before plugging it into an electrical outlet.

2. DANGER-VOLTAGE

- 1 Do not plug the charger base into an outlet with a voltage other than specified on the bottom of the charger.
- 2 This product is designed to operate within a range of 100-240 volts.
- 3 Voltage converters do not guarantee compatibility.
- 4 Do not force the plug into an outlet.
- 5 If the plug does not fit easily into the outlet, discontinue using.

3. DANGER-WARNINGS

To prevent the risk of burns, electrocution, fire, or injury to persons:

- 1 This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given adequate supervision or instruction concerning use of the product by an adult responsible for their safety. Children should be supervised to ensure they do not play with the product.
- 2 Use this product only for its intended use as described in this manual.
- 3 This product is designed to clean your teeth, gums and tongue only.
Do not use any attachments or brush heads other than those recommended by the manufacturer.

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This product has passed strict testing and inspection. Aquapick Co., Ltd. guarantees quality for 12 months from date of purchase.

We will repair defective products that occur under normal use within the warranty period free of charge, and replace the product with a new one if it is difficult to repair.

Precautions

Even during the warranty period, the following cases will be repaired at your expense:

- 1) The product is damaged due to improper use, mishandling or negligence after purchase.
- 2) The product is altered, modified or repaired by a party not authorized by Aquapick Co., Ltd.
- 3) The product is damaged as a result of connection to irregular voltage sources.

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